

LANAFORM

Limpia

BLACKHEAD SUCTION DEVICE



LANAFORM  
*Limpia*



## Complete skin care at home

The Limpia blackhead vacuum remover efficiently **removes pimples and blackheads**.

It efficiently fights skin impurities and promotes a healthy appearance of face skin.

Enjoy a personalised skin care thanks to its 3 intensities and extended battery life of 65 minutes.\*

\*Battery life measured at maximum intensity.








## Do you struggle to resist the temptation to squeeze your pimples?

Unlike squeezing pimples by hand, Limpia eliminates blackheads and impurities without risk of scarring or infection thanks to its **suction technology**.

It gets rid of imperfections faster and tightens your skin.

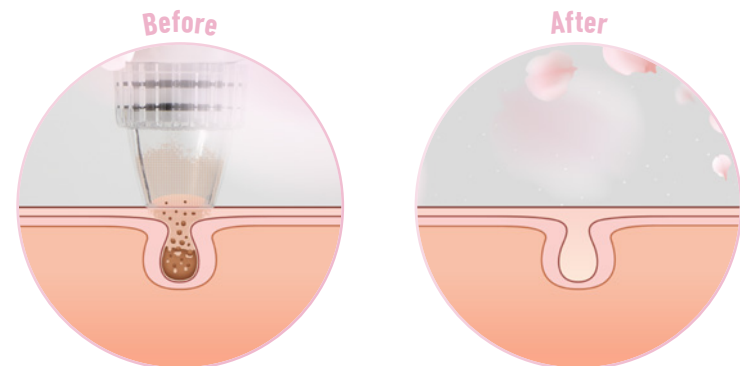

# Complete skin care at home

With five accessories, you will benefit from a complete skin care regime at home.



## Three different sized heads to remove impurities

Nose



The small head is perfect for hard-to-reach areas like the nose.

Cheeks & Forehead



The medium head provides an ideal suction surface for the cheeks and forehead.

Oily Skin



The large head is suitable for oily skin and removes excess sebum.

Anti-ageing



The elliptical head provides an anti-ageing treatment that smooths out wrinkles and promotes skin elasticity.


Dermabrasion



The Dermabrasion head gently exfoliates the face and removes dead skin.

## How to use

We recommend limiting the treatment to five minutes per use, with one treatment per week for combination and oily skin and one treatment every two weeks for normal and sensitive skin.

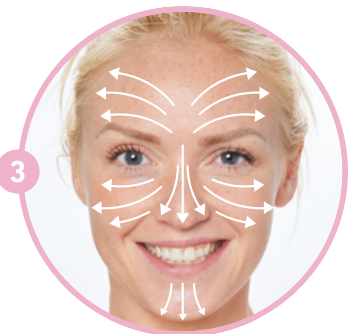

- 1 After cleansing your face, open your pores with a warm compress or facial sauna.
- 2 Select the tip of your choice and turn on the Limpia. If necessary, adjust the intensity before using it on your face.
- 3 Carry out your treatment.
- 4 Once the treatment is finished, cleanse your face thoroughly.
- 5 Apply your moisturiser or a mask to soothe your skin.

**Attention:** do not use the Limpia for more than two seconds in the same place and do not pass over an area you've already treated.

# LANAFORM *Limpia*

**LANAFORM**

## Features

- Deep pore cleansing.
- Five interchangeable tips.
- Three levels of intensity.
- Includes USB charging cord.



Master Cardboard Dimensions  
57.8 × 35.0 × 44.0 cm

Master Cardboard Weight  
17.000 kg

Colour Boxes / Master Cardboard  
50

Colour Box Dimensions  
21.0 × 10.7 × 5.8 cm

Colour Box Weight  
0.230 kg

Voltage  
5V





**LANAFORM**

Home Health Experience

Rue de la Légende 55  
4141 Louveigné, Belgium

Tel. : +32 (0)4 360 92 91

Fax : +32 (0)4 360 97 23

[www.lanaform.com](http://www.lanaform.com)  
[info@lanaform.com](mailto:info@lanaform.com)

